


# BodyZone<sup>™</sup>



## Advanced Pressure Redistribution Mattress

Advanced, innovative engineering and a revolutionary design provide excellent pressure redistribution from head to heel. The new BodyZone<sup>™</sup> patented technology is designed to provide optimal temperature and moisture management. The combination of our SMT<sup>™</sup> and VPF<sup>™</sup> help reduce friction and shear, while offering a higher level of envelopment to increase patient comfort.

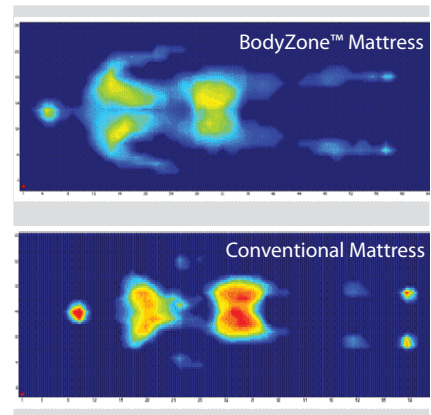

# BodyZone™ – Put your patients in the zone of full body comfort

## Prevention Mattress

- VPF™ TECHNOLOGY AND SMT™ GIVE OPTIMAL ENVELOPMENT AND UNSURPASSED PRESSURE REDISTRIBUTION AND PROVIDE SUPPORT FOR PATIENT WEIGHTS UP TO 500LBS
- SMT™ SIGNIFICANTLY REDUCES SHEAR AND FRICTION
- OUR OPEN CELL TECHNOLOGY PROVIDES ADVANCED TEMPERATURE AND MOISTURE MANAGEMENT WHILE BEING ENTIRELY ANTI-MICROBIAL




- EVERY BODYZONE MATTRESS IS MANUFACTURED USING AN ENVIRONMENTALLY FRIENDLY PROCESS (VPF™) AND USES RENEWABLE MATERIALS
- THE 360° STRETCH COVER\* MOVES IN EVERY DIRECTION AND IS IMPERVIOUS TO FLUIDS
- OUR MATTRESSES ARE AVAILABLE IN ALL SIZES AND INCORPORATE VPF™ TECHNOLOGY AND SMT™

### VARIABLE PRESSURE FOAMING PROCESS

- VPF™ is a technology used to create different levels of support with high envelopment and maximum comfort. This leads to better height and support retention over the life of the mattress.

### SURFACE MODIFICATION TECHNOLOGY

- SMT™ is the alterations in pattern, size of pattern, spacing between patterns and depth of cut to protect the most vulnerable areas of the body. This provides a zoned-cushioning layer, enhances air circulation and maximized comfort.

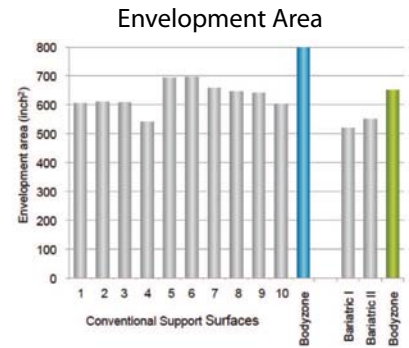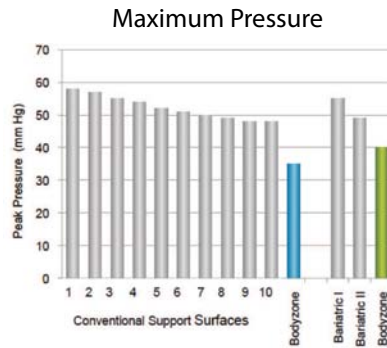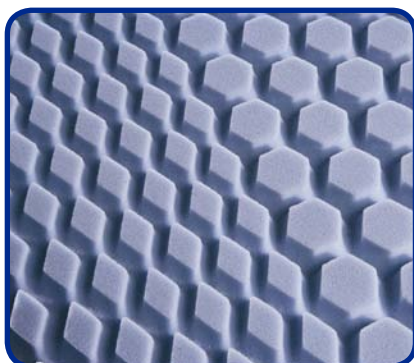

- When comparing the BodyZone™ mattress to a widely-used, regular standard polyurethane static mattress, the reduction of interface pressure was 49.2% less on the BodyZone™

- In lab testing vs. leading mattresses, BodyZone™ consistently yielded the highest area of envelopment (distribution of pressure across the mattress).



## Bodyzone Features & Benefits

### A PORTFOLIO OF PRESSURE REDISTRIBUTION MATTRESSES

	BodyZone™ 2000 New Standard	BodyZone™ 3000 Premium	BodyZone™ 5000 Premium Plus
360° Cover*	●	●	●
VPF™	●	●	●
SMT <sup>SM</sup>	●	●	●
Natural Ingredients	●	●	●
Open Cell Visco		●	●
Microspring <sup>SM</sup>			●

Microspring, an additional top layer of SMT foam, is added to further enhance patient comfort and reduce friction and shear.

Linet mattresses meet or exceed the requirements of 16 CFR 1632 and 16 CFR 1633.

I - Figures taken from a clinical abstract entitled »Pressure redistribution and envelopment with BodyZone™ mattresses: results of a series of tests.«

II - BODYZONE, SMT and VPF are registered trademarks of Foamex Innovations Operating Company.